
LINOS Microbench, Bench and Rail Systems

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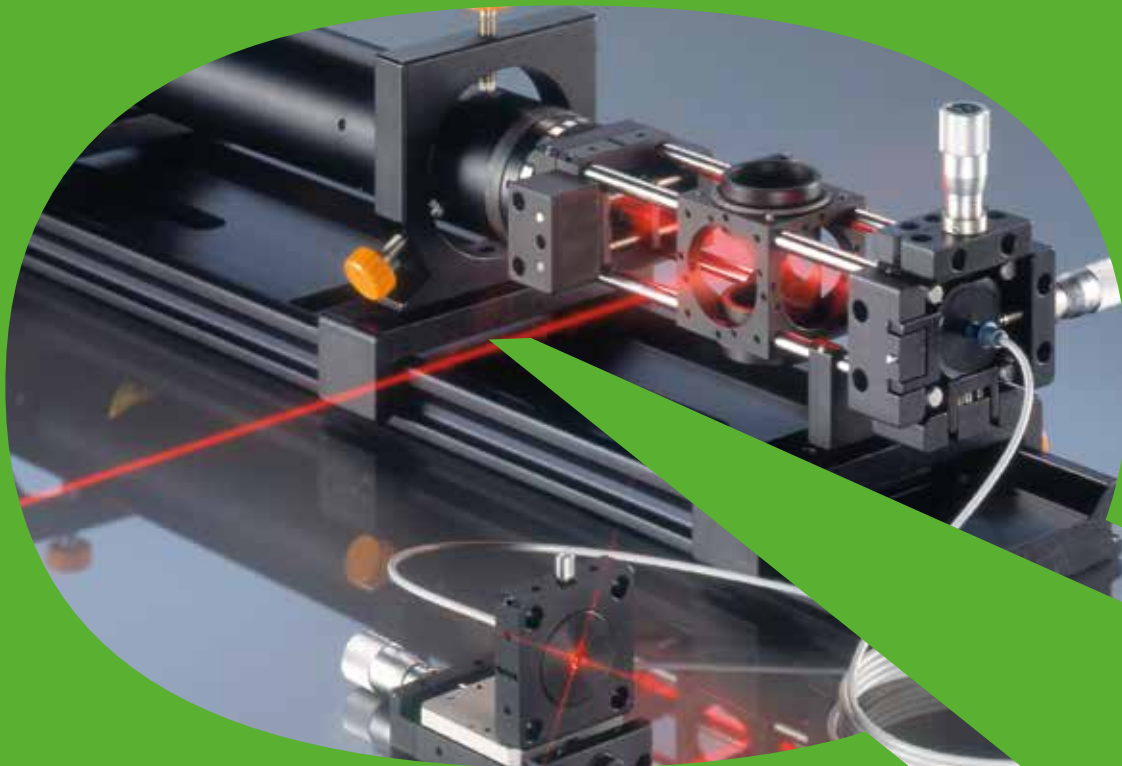NEW Edition 2018

About this Brochure

The LINOS Microbench™ system has been providing optical solutions for over 50 years. On this occasion, we are particularly pleased to present you a new edition of our Application brochure for LINOS bench systems.

Our primary focus is to provide you with a clear overview of the various systems as well as their compatibility with one another. Numerous photographs illustrate the variety of applications and setups, providing useful solutions for your applications.

Furthermore, this edition also contains our latest product offerings. For the first time we introduce our series of Q-Sets, compact optical subsystems built from standard catalog components. In addition, the brochure highlights the inter-compatibility of the LINOS Microbench™ with Qioptiq light sources and laser systems as well as Excelitas detectors.

Last but not least, you will find an application example in which the Microbench™ system enabled the rapid development of a complex optical measurement device - a flow cytometer prototype.

For further information please refer to the LINOS catalog as well as to our webshop, the Qioptiq Q-Shop: www.qioptiq-shop.com. There you can conveniently place orders online, view the latest products and receive the most up-to-date information.

Enjoy reading this brochure and don't hesitate to contact the Qioptiq catalog team for further technical advice at sales.qioptiq@excelitas.com or by phone at +49 551 6935-0.

1877

Rodenstock founded

1898

Spindler & Hoyer founded

1966

Pilkington PE Ltd. founded, which later becomes THALES Optics

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Content

| | |
|--|----|
| <i>The Microbench Concept</i> | 4 |
| <i>Microbench Manufacturing Line</i> | 7 |
| LINOS Bench Systems | |
| <i>First of its Kind: The Microbench</i> | 9 |
| <i>Microbench's Little Sister: Nanobench</i> | 9 |
| <i>Well-rounded Product Line: Tube System C</i> | 10 |
| <i>Getting on Track: Flat Rail Systems FLS</i> | 10 |
| <i>Basic Components of the different Bench Systems</i> | 11 |
| Frequently Asked Questions | |
| <i>How Do I Begin With My Experiment?</i> | 16 |
| <i>How Do I Integrate Spherical Optics?</i> | 18 |
| <i>How Do I Integrate Plano Optics?</i> | 19 |
| <i>How Do I Apply Adhesives?</i> | 20 |
| <i>How Do I Integrate Optics into Existing Setups?</i> | 21 |
| <i>How Do I Adjust Optics along the X and Y Axis?</i> | 22 |
| <i>What Options Are Available for Rotation of Optics?</i> | 23 |
| <i>How Do I Adjust Optics along the Z Axis?</i> | 24 |
| <i>How Can I Realize Different Angles in My Setup?</i> | 26 |
| <i>What Adapters are Available?</i> | 27 |
| <i>How Do I Mount My Assembly on e.g. Optical Tables?</i> | 28 |
| <i>What Light Sources Are Available?</i> | 30 |
| <i>How Can I Combine Systems?</i> | 32 |
| <i>How Do I Mount Bench Systems in My Instrument?</i> | 33 |
| <i>How do I connect Photon Counters to the Microbench?</i> | 34 |
| <i>How Do I Integrate Motorized Positioners?</i> | 36 |
| Q-Sets | |
| <i>Q-Sets – How to Save Time in the Lab</i> | 37 |
| <i>Q-Set Vertical Adjustment</i> | 37 |
| <i>Q-Set Variable Attenuator / Beamsplitter</i> | 38 |
| <i>Q-Set Spatial Filter</i> | 38 |
| <i>Q-Set Filter Cascade</i> | 39 |
| <i>Q-Set Fiber Coupler</i> | 39 |
| <i>Q-Set Spectroscopy Liquids</i> | 39 |
| <i>Q-Set Spectroscopy Plates</i> | 40 |
| <i>Q-Set Retroreflector</i> | 40 |
| <i>Q-Set Prism Stage</i> | 41 |
| Applications | |
| <i>Creating a Spatial Filter System</i> | 42 |
| <i>Creating an Expansion System</i> | 43 |
| <i>Creating a Fiber Coupler</i> | 44 |
| Prototype Setup | |
| <i>Setup of a Modular Flow Cytometer</i> | 45 |

2000

Rodenstock
Präzisionsoptik
acquired
by LINOS

2001

AVIMO Group
acquired
by THALES

2005

Qioptiq
founded as
THALES sells
High Tech
Optics Group

2006 / 2007

Qioptiq acquires
LINOS and Point Source
as "members of the
Qioptiq group"

2010

The new Qioptiq
consolidates all
group members
under one brand

2013

Qioptiq is acquired by
Excelitas Technologies

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The Microbench Concept

LINOS (the Catalog line of Qioptiq) – formerly Spindler & Hoyer – and the Microbench™ system are bound together inseparably, and have been so for more than 50 years now. The long-established Spindler & Hoyer company began production of the compact 4-rod Microbench™ system, the innovative successor to the triangular rail system, as early as 1968.

Until well into the 1960s, optical test assemblies were typically set up on triangular profile, or "Zeiss profile", optical benches. These were heavy and relatively large constructions that permitted only linear, level optical arrangements. Those systems that did allow spatial expansion were not sufficiently torsion-resistant for metrological applications and thus were only suitable for demonstration purposes. The dimensions of the accessories for these earlier model benches made it impossible to achieve the compactness often sought-after in an optical setup. Furthermore, the insertion of microscope assemblies such as polarization, interference or measuring microscopes was exceedingly complex with conventional optical benches. The logical step was the development of a novel micro-optical bench.

The distinctive feature of the new system was its high stability and precision. Since its introduction, optical and mechanical components of the Microbench™ system have become indispensable in optical laboratory experiments and are built into a large number of optomechanical devices.


Instruments made up of Microbench™ components are also ideal for demonstrations and practical training: the exposed design makes it not easy only to observe the entire system but also to alter or expand its most basic workings by simply adding components, such as grids, prisms, mirrors or other accessories.

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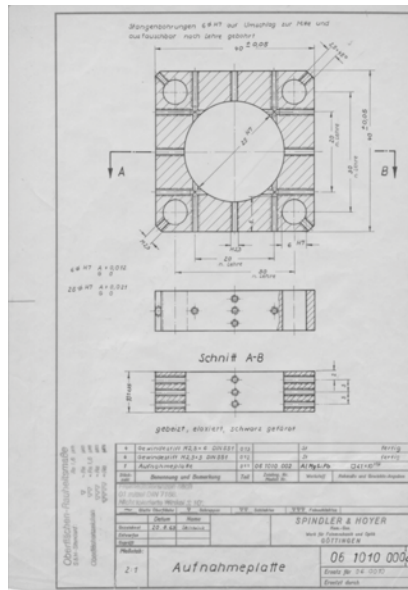In the following years the Microbench™ system has been continually expanded. Numerous mounted optics, specially designed for the Microbench™, have been developed. The special  icon has been introduced for use in our catalogs to make it easier for you to find mounted optics to add to your Microbench™. Now in use both in our main catalog and the on-line shop, this icon shows you at a glance which optics are available in Microbench™ compatible mounts.

For more than 50 years the Microbench™ system has maintained compatibility of all elements. Thanks to our extremely strict production tolerances, we can guarantee that components made today can still be used with elements manufactured in the 1970s. In our further development of the Microbench™ system, we will continue to ensure that new products are compatible with earlier-model elements.

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Your advantages

- Decades of experience
- High flexibility thanks to a wide range of components
- High quality ensured by German manufacturing standards

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Microbench Manufacturing Line

Top-quality materials

- RoHS compliant materials
- Special alloys

meet the highest standards

Optimized manufacturing

- Parts manufactured in a single work cycle
- Highest precision to meet our stringent tolerances

for pinpoint accuracy

100% Quality control

- Final inspection of all tolerances
- All components checked for accuracy of fit
- Compatibility of all Microbench components assured

for consistent quality – for more than 50 years!

Surface finishing with precision anodizing

- Resistance to wear and tear
- Highest abrasion resistance ensures stringent manufacturing tolerances
- Matte black anodizing prevents reflections


for best surface quality

Materials made from highest-quality aluminum alloy, for example: bulk material (round) for mounting plates

9-axis CNC milling machine handles all production stages, from aluminum bulk material (round) to finished mounting plate with all drill-holes and threads

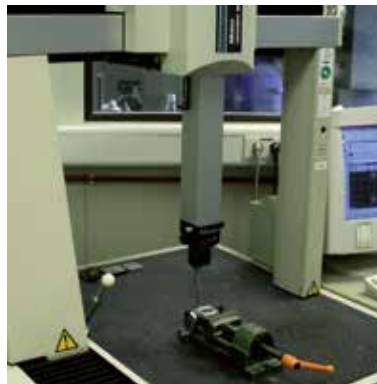Final inspection of mounting plate tolerances at a coordinate measuring station

Microbench cubes after surface finishing in the anodizing bath

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Rail System FLS 40

LINOS Bench Systems

*LINOS bench and rail systems for experimental set-ups and measuring / inspection equipment offer high precision, excellent mechanical stability and flexibility in use. Our **Microbench** system has stood the test of time for more than five decades, and we are constantly expanding the system to make sure all your requirements are met. At the same time, we continue to specialize in trend-setting developments.*

*Several years ago, in response to the growing demand for smaller systems, Qioptiq introduced the **LINOS Nanobench** system. Compared to the 40-mm edge length of the Microbench system, the Nanobench edge length is just 25 mm.*

*Another product line is our modular, closed component **Tube System C**. It enables fast set-ups of enclosed optical systems. The tubes have outer diameters of 30 mm and 35 mm. The system features a wide variety of lens mounts, diaphragms and adjustment units. The universal 1-32 UN mounting unit (C-Mount) enables direct connection to CCD cameras, lasers and microscopes.*

*When searching for mechanical systems for micro-optical applications, look no farther than to our **Rail Systems FLS**, including the FLS 40 series specially developed for the LINOS Microbench. This flat rail system can be used with individual components, such as mounted optics or positioners. It can also be used to fasten complex Microbench structures to optical tables.*

No matter which products you choose from our modular bench and rail systems, excellent compatibility of the systems with one another is guaranteed. And many other of our products, from mounted optics to lasers, can be directly integrated as well. This gives you the advantage of fast implementation, whether for experimental set-ups or permanent installation of complex measuring equipment. Whatever your photonic application is, you can find the optimum system or combination of systems for your particular specifications and spatial conditions within the extensive range of our bench and rail systems.

☒ **Microbench and Nanobench:**
optomechanical bench systems

☒ **Tube System C:**
modular and closed component system

☒ **Rail System FLS:**
flat rails and carriers for high-stability structures

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- Ideal system for optics with \varnothing 18 to 31.5 mm
- Modular design ensures greatest flexibility
- Extensive range of optics in Microbench mounts


First of its Kind: The Microbench

Microbench is a compact component system for numerous areas of application – from lighting and imaging structures to optical experimental setups, measurement / inspection equipment and interferometric applications, to name just a few.

The outstanding characteristics of the LINOS Microbench products include their high mechanical stability - thanks to the use of highest quality materials - , as well as their high axial precision due to narrow production tolerances. The Microbench is easy to handle, even for first-time users.

A high degree of flexibility is ensured not only through the wide variety of components available for LINOS Microbench systems, but also by the range of adjustment elements and the system's modular design.

The basic principle lies in the combination of standardized mounting plates for integrating optics of all types, connected by stainless steel rods. Both linear and spatial three-dimensional setups can be implemented quickly and easily.

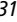We offer mounted optics developed specially for the Microbench with diameters from 18 mm to 31.5 mm. The special  icon, seen in both our main catalog and the online shop, shows you at a glance which optics are available for your Microbench system.

Microbench's Little Sister: Nanobench

- Ideal system for optics with \varnothing 3 to 12.7 mm
- Modular design ensures greatest flexibility
- Extensive range of optics in Nanobench mounts

The Nanobench is an extremely compact construction system for optics with diameters of up to 12.7 mm, and features exceptionally high precision. As with the Microbench system, mounted optics are fastened in mounting plates, which are connected by four rods to ensure extremely precise adherence to the optical axis.

A variety of base plates make it easy to fasten Nanobench components directly onto optical tables, mounting plates and support rods. With its compact dimensions and flexibility in three-dimensional structures, the Nanobench is ideal for applications in micro-optics and fiber optics. Here, too, Qioptiq offers a broad selection of mounted optics that can be directly attached to the Nanobench. The  icon helps to make it easy to find Nanobench-compatible products and accessories when ordering

Adapters are available for fitting Microbench and Tube Mounting System C elements to your Nanobench system (please look on page 27).

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- For fast realization of your laboratory instruments, samples, prototypes and small-batch productions, including lenses and expander systems
- Modular component system
- For unmounted optics with \varnothing 6 to 31.5 mm
- Compatible with C-Mounts, microscopes and CCD camera threads

Well-rounded Product Line: Tube System C

The Tube System C is a modular component construction system for fast and trouble-free set-up of laboratory instruments, prototypes and small-batch production. Tubes are available in diameters of 30 mm and 35 mm. Threaded sockets and threaded rings enable the connection of a wide range of optics, with diameters from 6 mm to 31.5 mm

With the Tube System C, creating your own lens and expander systems is a breeze. Adapter rings make this system compatible with the connecting threads on C-Mounts, microscopes and CCD cameras. Fine adjusters and focusing mounts let you make small changes for highly precise focusing. The Tube System C and Microbench are excellent in combination; for example, to create a light-proof and dust-tight beam path.

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- FLS 40 compatible with LINOS Microbench
- For mounted optics with \varnothing 25/30/35 mm
- Direct installation on breadboards, optical tables, etc.
- Great flexibility and compatibility

Getting on Track: Flat Rail Systems FLS

Our flat rail system FLS includes the FLS 40 and FLS 65 types.

With the FLS 40, you can implement linear and areal structures and fasten mounted optics with diameters of 25, 30 or 35 mm, in conjunction with mounting plates or cubes from the Microbench series fixed onto a carrier. Because FLS 40 is fully compatible with the LINOS Microbench system, you can add basic components and other special elements such as mirror and prism adjusters, positioning units and fiber couplers. The FLS 40 has an optical axis height of 40 mm.

Flat rail system FLS 65 as well as the Profile system X95 are suitable for set-ups in larger dimensions. Those systems can also be combined with the Microbench with the help of adapters.

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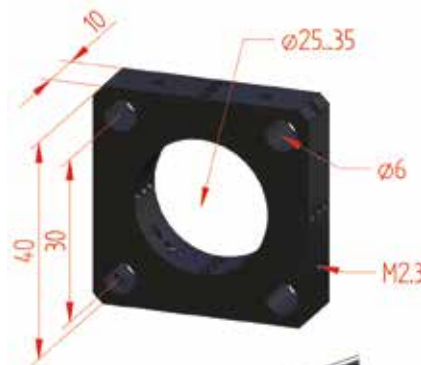
Basic Components: Microbench

The Microbench system is based on a small number of elements, such as mounting plates, mounts and cubes, in which optics are installed. These mounts are connected by rods to form extremely compact and sturdy units.

The bore holes for the rods on all **mounting plates** have identical outer diameters and spacing. The edge length is 40 mm, and thickness is 10 mm. Where not indicated otherwise, the same dimensions apply for other components in the series, such as holders and special mounting plates.

- Lead-free aluminum alloy
- Matte black anodized surfaces
- Inner diameter: 25, 30 or 35 mm with H7 tolerances
- Four bore holes for mounting the connecting rods. The tolerances for the rod bore holes on the mounting plates, together with the tolerances for the rods, guarantee a system with virtually zero play.

Stangen / Rods

Stellringe G061220000 (10 St.)
Set of collars G061220000 (10 pcs.)

The **rods** are made of corrosion-resistant stainless steel and are specially designed for connecting individual Microbench components. They are induction-hardened to prevent marks made by adjustment screws. To ensure that there is no play in the guidance of mounting plates and other components, the rods are lapped to fit.

- Hardened stainless steel with polished surfaces
- Lapped to fit
- Outer diameter: 6 mm
- Lengths: 20 mm to 450 mm
- Compatible collars as holders or stops; e.g., for Z-adjustment made of individual components

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Basiskomponenten / Basic components

Basic Components: Microbench

Cubes are multifunctional elements for steering light beams and for integrating various types of optics in holders. Cubes can be combined with rods and mounting plates to create compact, three-dimensional structures. They are available with mounting diameters of 25 mm and 30 mm to fit holders in different sizes.

- Multifunctional elements for 3D structures
- Ideal for beam deflectors
- Inner diameter of bore holes: 25 or 30 mm
- Two of the cube sides have bore holes for rods, four sides have threads for mounting plates with tapped holes

For a large assortment of mounted optics made specially for the Microbench, please see the LINOS catalog or www.qioptic-shop.com (look for the Microbench icon: )

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Mounting accessories


Microbench screws and threaded pins have M2.3 threads and are essential accessories for all mounting plates, holders, and cubes.

Please find detailed information regarding the complete Microbench accessories here:

- LINOS Catalog, Edition λ , Page 75
- Q-Shop <https://www.qioptic-shop.com/en/Optomechanics/LINOS-Microbench/Tools-Screws/>

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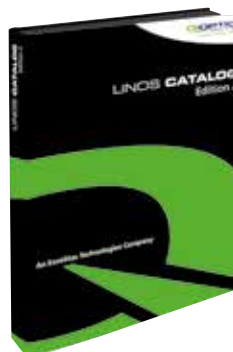LINOS Catalog, Edition λ

| Order No. | Product | Qty. |
|---|------------------------------|----------|
| G06 1011 000 | Set of threaded pins, M2.3x3 | 150 pcs. |
| G06 1012 000 | Set of threaded pins, M2.3x6 | 150 pcs. |
| G06 1090 000 | Screwdriver | 1 pcs. |
| Alternative: For fastening optics in mounts without tools | | |
| G06 1238 000 | Set of thumbscrews, M2.3x4 | 20 pcs. |
| G06 1239 000 | Set of thumbscrews, M2.3x9 | 20 pcs. |

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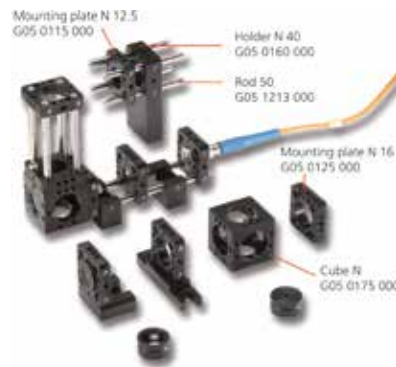
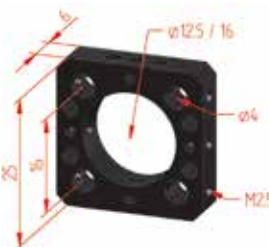
Zubehör / Accessories

Basic Components: Nanobench

The functional principle of the Nanobench is identical to that of the Microbench system. The dimensions and tolerances, however, are on a smaller scale, compared to the Microbench's size. A wide range of mounting plates, holders, rods, and cubes are available for creating sturdy and extremely compact systems.

The bore holes for the rods on all Nanobench mounting plates have identical outer diameters and spacing. The edge length is 25 mm, and thickness is 6 mm. Where not indicated otherwise, the same dimensions apply for other components in the series.

Mounting plates

- Inner diameter: 12.5, or 16 mm with H7 tolerances

Rods

- Outer diameter: 4 mm
- Lengths: 12 mm to 200 mm

Cube

- Inner diameter of bore holes: 16 mm
- To mount beamsplitters, mirrors and prisms with mirror holder 45° or prism holder
- Connection from cubes to rods with mounting plates N 16
- Set of screws (G05 1522 000) are included

Please find detailed information regarding accessories:

- LINOS Catalog, Edition λ , Page 99/100
- Q-Shop <http://www.qioptiq-shop.com/en/Optomechanics/LINOS-Nanobench/Tools-Screws/>

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Tuben / Tubes

Basic Components: Tube System C

Key components of the tube mounting system are the Tubes C, which are available in two diameters and various lengths. All tubes are made of matt black anodized aluminum.

Tubes C30 have an outer diameter of 30 mm and a free inner diameter of 22.4 mm. Each features C-Mount connection threads inside and outside. Tubes C35 have an outer diameter of 35 mm and a free and are available in lengths from 20 mm to 80 mm. The C35 tubes feature an M32.5x0.75 inner thread for mounting optics.

- Light-tight and dust-tight system for use in research, development and industry
- For fast and economical set-up of customized optical systems (e.g. prototypes and small-batch production)
- Universal mounting unit enables direct connection to CCD cameras, lasers, microscopes
- Compatible with Microbench

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Beam deflector Cube C lets you connect mirror and prism accessories. These have 30 mm bore holes for prism accessories and three (or five) bore holes with C-Mount inner threads for connecting the tubes

Optikmontage / Optics mounting

Unmounted optics can be installed directly in the Tube system C and secured using threaded rings.

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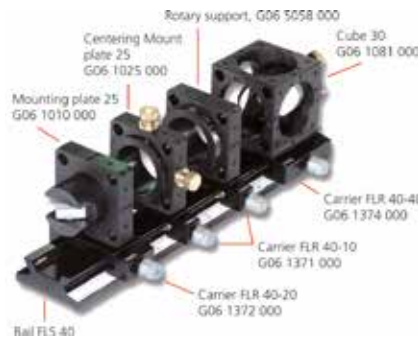
FLS 40

Flat rail System FLS 40 is compatible with our Microbench systems with the corresponding dimensions and beam heights. The system features excellent straightness and evenness with low torsion.

The slotted holes in the rails enable mounting on breadboards and optical tables independent of existing hole patterns.

The carriers have a spring-loaded lever. Just loosen the knurled screw, and the spring tension makes it easy to move the carriers smoothly to the desired position. Then tighten the knurled screw again to affix the carrier in position.

Precision milling of bearing surfaces together with abrasion-resistant anodizing ensure that components slide smoothly and evenly on the rail.

- Precise bearing surfaces
- Abrasion-resistant anodizing
- Lengths: 40 mm to 1000 mm
- Scale
- Slotted holes on the rail enable table mounting independent of existing hole patterns
- Compatible with Microbench system
- Carrier with spring-loaded clamp
- Optical axis height of 40 mm

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FLS 65

The flat, extruded rail features excellent straightness properties and low torsion. Bearing surfaces are precision processed for smooth movement of the carrier.

The rail features a matte black anodized, abrasion-resistant surface. Slotted holes make it easy to fasten the rail to optical tables independent of existing hole patterns, using either M6 or 1/4" screws. The FLS 65 rail is available in lengths from 200 mm to 1500 mm.

The carriers FLR 65 are available in three widths and feature a universal bore hole pattern for attaching optical mounts. The spring-loaded clamp ensures smooth movement of the carrier along the rail, and a knurled screw fixes the carrier in position. Just like the rails, the carriers feature matte black, abrasion-resistant anodizing. Together with the rail, the surface of the carrier is 25 mm above the base.

The optical axis height is 65 mm above the working surface when standard optic holders for the FLS 65 system are used (please see the mounting example).

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Simulation of the beam path of a basic alignment telescope with Qioptiq WinLens™

For detailed information and a free version of the WinLens™ 3D software have a look at: www.qioptiq-shop.com/Optik-Software/Winlens-Optical-Design-Software/

We offer WinLens Webinars in English with our Specialist. Please have a look regarding available training dates.

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
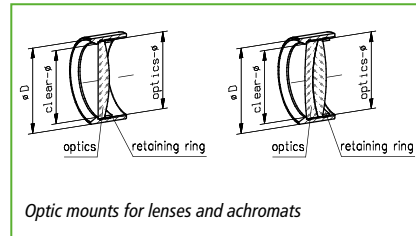
Optic mounts for lenses and achromats

Required optics in different mounts depending on the outer diameter of the lens

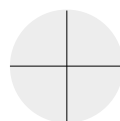Reticle crosshair, unmounted G391130000

How Do I Begin My Experiment?

The following pages serve as a checklist for a structured planning and design of your setup, to help you reach an optimum, economical solution for your requirements. This example describes the construction of a basic alignment telescope.

How do I design my optical assembly?

We recommend using an optical design software, such as WinLens™ 3D from Qioptiq, to create a simulation. This can be a tremendous help in preparing your optical assembly. This software lets you simulate the desired beam path and analyzes the results directly. You are able to optimize the setup as needed before finalizing your design.

Furthermore, WinLens™ 3D comes with an integrated database that contains all of our standard optics. Using standard optics helps to keep costs to a minimum, and mounted optics for direct integration in our construction systems are also available.

What optics do I need and which construction system is the best?

Your optic diameters and areas of application will play a role in determining the optimum construction system for your requirements.

In our example, the Microbench system was chosen. In this example, we have decided not to enclose the beam path. As an alternative, we can recommend the use of the Tube System C for a closed construction.


The following optics are required for this project: A plano-convex lens $f=80/22.4$ (G312317000), a reticle 19 mm (G391130000), a second plano-convex lens $f=40/22.4$ (G312314000) and an achromat $f=60/25.4$ (G322306322).

For mounting in the Microbench mounting plates, there are suitable mechanical mounts for the optics and lenses: Mount CL 19 (G063661000), Mount CL 25.4 (G063663000) and 2x Mount CL 22.4 (G063651000).

The combination of this choice of optics for the alignment telescope creates an inverse image with a 3x magnification. The distance setting (focus) is adjustable in the range of 300 mm to ∞ , with the built in line reticle plate as a useful alignment aid when measuring distance deviations with respect to a reference line.

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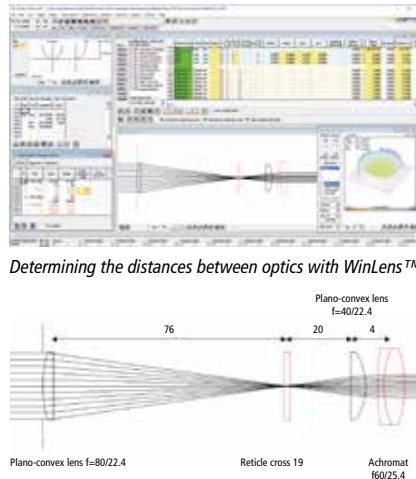How do I position the optics? And what about the adjustment?

After the beam path has been laid out and the optics have been selected, the distance between each of the optics must be determined. You can either calculate the distances yourself or use the Qioptiq WinLens™ program for this purpose. In this example, WinLens™ offers the advantage that it determines the required adjustment range on the Z-axis using its online optimization feature. This is important for your choice of adjustment components.

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
The mechanical components

Which mechanical components do I need?

In addition to the above-mentioned mounts for lenses and optics, further components are now required for a suitable opto-mechanical framework for setting up the lens system. The following solution is based on the Microbench system.

Mounting plates 25 (3x G061010000) and a mounting plate 30 (G061042000) are used for fixing the mounted optics. They are aligned on rods (2x 150 mm, G061210000) and can be fixed directly on an optical table via standard rod holder F (2x G061210000).

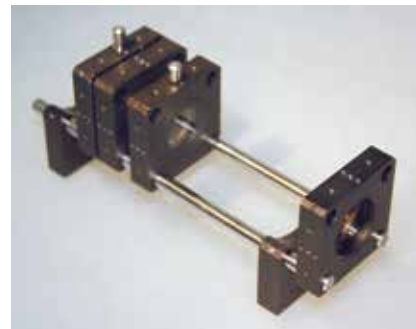
The finished setup of a basic alignment telescope

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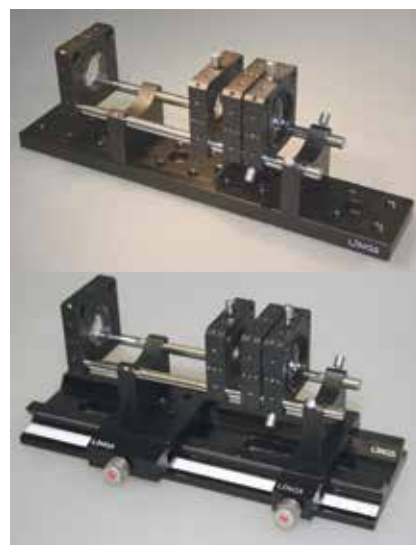Above: The setup mounted on a Mounting plate 220x40
Below: The setup mounted on a Flat rail FLS 40

How should my structure be secured?

We offer you a broad range of equipment and methods for fastening your construction to a work surface. You can fasten it to a mounting plate using screws, mount it on an optical table or breadboard, or integrate it into an existing instrument of your own.

For answers to these questions, please see pages 28-29, respectively page 33.

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For a large assortment of mounted optics made for the Microbench and Nanobench, please see the LINOS catalog or www.qioptiq-shop.com (please have a look for the icons  and )

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Unmounted optics ø 31.5 fixed with retaining rings M32 5x0.75, G03 6051 014
Unmounted optics ø 22.4 fixed with retaining rings M23 2x0.75, G06 3666 000

How Do I Integrate Spherical Optics?

Mounted spherical optics

The easiest way to integrate spherical optics in the Microbench system is by using mounted optics from our extensive product portfolio. Depending on the particular optics, the mounts have an outer diameter of 25, 30 or 35 mm and are ready for installation in the compatible mounting plates, holders and adjustment elements. The mounts are fastened by using threaded pins.

Mounted spherical optics

For the Nanobench system, a wide range of mounted optics is available. The outer diameter of the 16 mm mounts corresponds to the mounting diameter of the Nanobench mounting plates, holders, cubes and adjustment elements.

Lens mounts

Of course, lens mounts for the Microbench and Nanobench system can be purchased separately. We offer Nanobench mounts for optics diameter from 3 up to 12.7 mm and Microbench mounts for 18 up to 25.4 mm.

Unmounted spherical optics

The Microbench system also offers the option of integrating unmounted optics. For unmounted spherical optics with diameters of 25 mm, 31.5 mm, 0.5" and 1" we offer mounting plates with M35 threads and compatible threaded retaining rings.

Optics for the tube system

The tube system is constructed for mounting unmounted lenses and achromats with diameters of 22.4 mm (in Tube C30) or 31.5 mm (in Tube C35). The optics are fixed in the desired position with threaded rings.

For mounting smaller optics with an outer diameter of 6 to 19 mm, you can use threaded mounts. Regular or achromatic lenses can be fixed in the compatible threaded mount which has outer threads for installing and positioning within the tube. The high precision of the overall system ensures that lenses are centered along the optical axis.

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☑ With adjustment options

If you want to mount prisms or beamsplitter cubes in Microbench cubes, we recommend Adjustable prism insert 30. For oval mirrors and beamsplitter plates, we recommend Adjustable mirror insert 30. To mount prisms, mirrors and cubes in mounting plates, use Tilting platform 25. Prism table 25 comes installed in a mounting plate and can be rotated.

Qioptiq also offers a fully-equipped Beam steering mirror assembly, including front-surface mirror and housing, specially made for adjustment of elliptical mirrors. It is also available with a UV mirror.

On all adjustable optics holders, the optics have to be bonded to the desired holders.

Please find straight tips regarding working with adhesives on page 20.

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
☑ Without adjustment options

The Clamp holder for prism enables fast and simple mounting of prisms and cubes with a 10 mm edge length. Prism table 25 is the ideal mount for edge lengths of up to 20 mm. On both components, the optics are secured by clamps.

The Optic support 25 is used for mirrors, which are bonded to the guide edge. These mirror mounts can be installed in Mounting plate 25 and in Cube 25.

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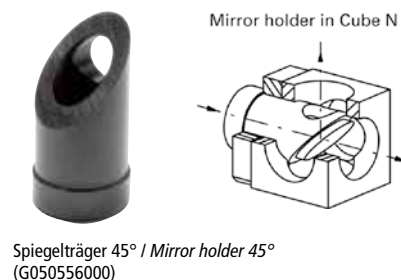☑ For elliptical mirrors of 12.5x18 mm, the Mirror holder 45° is the ideal mount, while the Mirror mount 45°/ 7x12 is best for beamsplitter plates with dimensions of 7x12 or a mirror diameter of 10 mm. Prisms and cubes that measure 10x10 mm are bonded to the prism holder.

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Applying adhesive in spots between optics and contact surface of the holder

Positioning paper strips as spacers on a mirror mount

Pressing down and repositioning a mirror on the fixed mirror mount

Beamsplitter cube illuminated by a UV light source

- To keep the need for cleaning to a minimum, wear latex gloves or lint-free cotton gloves when unpacking and setting up the equipment and make sure all work surfaces are clean.

How Do I Apply Adhesives?

Optics:

Use a bulb blower to remove dust and loose dirt. Put a lens-cleaning cloth or another soft cotton cloth (e.g. Opto-Wipes™) on the optics and use a pipette to apply several drops of alcohol to the cloth. Make sure the alcohol is chemically pure. Then drag the cloth slowly over the optics. To reduce fingerprints or similar impurities, proceed as described above or soak a cotton swab in alcohol and wipe it across the optics in an "S" shaped motion.

Optics holder:

The holder must also be cleaned and freed of dust and grease. To do this, wipe it with a soft cloth, using a small amount of a grease-cutting detergent (such as dishwashing liquid) if necessary. Wipe again afterwards with a fresh cloth to remove any residues.

Adhesives:

The type of adhesive you should use depends on the type of optics you wish to cement. For mirrors, use an elastic two-component adhesive for a uniform curing and tension-free hold. For components such as beamsplitter cubes or plates, a UV adhesive is ideal. As the name indicates, UV light is applied to cure the adhesive.

Strips of double-sided tape are commonly used in research and development laboratories. These are well-suited for quick experimental setups, particularly if the optic is removed from the carrier afterwards. They are not recommended, however, for long-term use; due to their elasticity and susceptibility to vibration, it is impossible to ensure that the position of the optics is maintained.

Tips for cementing optics:

Even if a liquid adhesive is used, the optics can be removed later if the adhesive is applied in discrete spots on the edge of the optics only. Affix the optic holder to your work surface with adhesive strips, making sure the glue surface is on the horizontal plane to prevent slippage until the adhesive is cured. To ensure an even layer of adhesive, lay two strips of paper as spacers on the optics carrier and press the optics until the adhesive reaches the paper. Once the adhesive has cured, you can remove the paper. This method also prevents the adhesive from leaking out between the optics and the holder when the parts are pressed together. The use of a cotton swab is recommended to prevent damage when aligning and pressing down on the optics.

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- Use magnetic mounting plates for quickly switching optics.

How Do I Integrate Optics into Existing Setups?

☑ In addition to standard mounting plates as basic components, Qioptiq offers snap-in (latchable) mounting plates and a magnetic version. They are ideal for adding optical components to an existing system, and make it possible to expand a system without having to remove the rods. With the mounting plates, for example, you can quickly replace optics such as lenses, filters, reticules, and others in the beam path. For unmounted optics with diameters of 1/2", 1", 25 mm and 31.5 mm, we offer the M35x0.5 snap-in mounting plate, in which the optics are mounted with compatible threaded rings.

By mounting the snap-in (latchable) plate in such a way that it can be swiveled, optics can be removed from the beam path without being removed from the setup. Collars are inserted to fix the position of the plate.

Magnetic mounting plates are easy to fasten between components on a four-rod system. The magnets are strong enough to attach the mount even upside-down which make them a comfortable accessory for your optical setup.

☒ For the Nanobench, there are also a number of options for retrofitting your assembly. Snap-in (latchable) and swiveling (rotatable) mounting plates are the perfect solution for adding mounted optical components to an existing system without having to remove the rods.

With these plates you can quickly switch between different optics in experimental setups. This type of mounting plate is also ideal if you want the option of removing optics from the beam path without removing them from the assembly completely. Simply install the mounting plate so that it can swivel. Collars are inserted to determine the position of the plate.

- Mounting Plate N, snap-in and swiveling versions
- For retroactive installation in 2-rod, 3-rod and 4-rod structures
- For mounted optics with \varnothing 16 mm

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How Do I Adjust Optics along the X and Y Axis?

Q For adjustment in two spatial directions, the X-Y Precision translation stage is ideal.

The Centering mounting plate 25 (also available as a holder) is equipped with easy-to-turn brass screws with a 0.25 mm pitch, permitting an adjustment range of ± 1 mm. The rotatable inner ring enables adjustment, for example, of crosshairs.

The most precise X-Y adjustment is provided by the fle.x plate with an adjustment range of ± 0.25 mm. It features a solid flexure and excellent longterm stability. The fle.x plate is best suited for coupling into a light guide, for example an optical fiber. The mounting diameter is 12.5 mm.

Q Economical yet precise adjustment is possible with the Adjustment plates for the Microbench and Nanobench systems. Two plates connected by springs are adjusted to one another and fixed with two screws. The Slip plate positioner M for the Microbench is available with mounting diameters of 25 mm and 30 mm.

N The X-Y Translation stage N is ideal for use in the Nanobench system. The mounting diameter is 16 mm for mounted Nanobench optics, and the mount has a clear aperture of 8 mm. The fine-pitch screws have a 0.25 mm pitch and an adjustment range of ± 0.5 mm.

■ The Spatial filter module C30x37.5 is specifically designed for positioning pinholes. The solid flexure enables a high-precision adjustment that is free of backlash. C-Mount connections on both sides and an outer diameter of 30 mm enables versatile use with the Tube mounting system, with the Microbench, and on the FLS 40 flat rail. Mounted pinholes of various diameters are available to fit into M16 inner threads.

You want to construct or document the setup in your CAD software?

On www.qioptiq-shop.com we offer CAD files in dxf or step format. You will find the files for the respective products under the tab Docs + Drawings

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What Options Are Available for Rotation of Optics?

☑ We offer two Rotary supports for mounting components that require rotational alignment, such as slits, reticles and polarizers. The first is a mounting plate with an integrated rotary mount, which features a scale with 5° intervals and a rotation range of 360°. The inner ring can be fixed with a locking screw.

For higher precision requirements we recommend the Rotary support 20°, with a micrometer screw for adjustment. The adjustment range is $\pm 10^\circ$ with a sensitivity of 10". Both rotary supports have mounting diameters of 25 mm.

☑ Motorized rotation mounts are specifically designed for fast rotation of construction elements with 25 mm diameter. You can choose from either servomotor or stepper motor models. In the servomotor model, the motor speed can be adjusted from 10 to 200 rpm. The maximum speed of the stepper motor model is 400 rpm, and the smallest step angle is 1.8°. The servomotor requires a supply voltage of 0-12 V for activation, while the stepper motor requires a controller. An optional reflected-light barrier is available separately.

☑ For the Nanobench system the Rotary support N is most suited. This component is designed to be clipped into Mounting Plate N 16 and is tensioned by a disc spring to ensure a rotation without backlash. A lateral bore hole with a diameter of 1.1 mm is provided for insertion of an optional lever to facilitate turning the support.

The Rotary support N can also be used as a quick-change adapter for optics in mounts with a diameter of 16 mm, as the mounted optics are simply clipped into place.

- Rotary Support for 360° rotation and Rotary Support 20° for more stringent requirements
- Motorized rotary mount with stepper motor or servomotor

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How Do I Adjust Optics along the Z Axis?

☐ You can put together a simple Z-Translation stage made up of a Holder for micrometer, a Micrometer screw, the Mounting plate to be adjusted, Compression springs (2 pcs.) and Collars. You can also use a different construction element instead of the Mounting plate.

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☐ The Z-Adjustment M is made up of a Fine adjustment screw, a Compression spring and a Collar for high-precision fine adjustment within a 5 mm range. It enables a backlash-free adjustment of optics, pinholes, fiber connectors, and more. Since the Z-Adjustment M is extremely compact it can be retroactively integrated into any structure.

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☐ For an extremely high precision adjustment along the Z axis, we recommend the Z-Fine Adjustment M. It enables direct mounting of components with a diameter of 25 mm. The micrometer screw with 0.25 mm pitch ensures backlash-free and high-precision adjustment along the Z axis, with a 1.5 mm range and a 1 µm sensitivity. Its integrated 1-32UN inner threads make it compatible with the Tube Mounting System C.

Because it can be accurately centered with virtually no tilt, this precision positioner is well suited, for example, for beam coupling in fibers.

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- Simple Z-Translation stage made of Holder with micrometer, Compression springs, Collars and Mounting plate
- Z-Adjustment M for precise fine adjustment with 5 mm adjustment range
- Z-Fine Adjustment M with a 1.5 mm range enables high-precision and backlash-free adjustment of optics, pinholes, fiber connectors, and more

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The Z-Adjuster N consists of a Fine adjustment screw, a Collar and a Compression spring. You can combine it with the Mounting Plate N to create a simple Z-Translation stage. Instead of a Nanobench Mounting plate, you can use the Holder N 40 or Adjustment plate N.

The Z-Adjustment N has a range of 5 mm and provides a backlash-free adjustment of optics, pinholes, fiber connectors and more. It is extremely compact and can be retroactively integrated into Nanobench structures.

Our Fine adjustment C30 and Focusing mount C35 let you alter the spacing between Tube System C elements along the Z-axis.

The Fine adjustment C30 is available with an adjustment range of 3 or 6 mm (with different thread pitches). The overall length of the Fine adjustment C30 changes with adjustment.

The Focusing mount C35 enables the linear adjustment of optics in threaded mounts or of directly installed optics with a diameter of 22.4 mm, without changing the overall length of the mount. If you wish to adjust the zoom lens in an objective from individual components, the Focusing mount C35 is the ideal component.

Both, the Fine adjustment C30 and the Focusing mount C35 feature C-Mount threads and can be fitted directly between the tubes.

- Fine adjustment C30 for Z-adjustment of optics with a range of 3 mm or 6 mm
- Focusing mount C35 for linear adjustment of optics in threaded mounts. The total length of the mount remains constant

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You want to construct or document the setup in your CAD software?

On www.qioptiq-shop.com we offer CAD files in dxf or step format. You will find the files for the respective products under the tab Docs + Drawings

How Can I Realize Different Angles in My Setup?

There are a number of options available for setting up angle connections at other than right angles. One is the Connector 80 V, which connects two mounting plates on one level. With this connector you can have one beam branch off diagonally with respect to the axis. The slotted holes enable the continuous adjustment of the angle and the mounting plates. The crossed struts help you achieve an optimum stability.

Connector 80 V is also ideal for fixing the angle of the swivel connector. The Swivel connector ensures a connection between two mounting plates with no play, and has a pivot range of 120°.

- Connector 80 V for creating variable distances and angled connections
- The Swivel connector has a swivel range of 120° with no play
- Swivel arm with scale for measurable swiveling of Microbench assemblies through $\pm 90^\circ$; direct installation of mounting plates or Rod holder F

Another option for setting up a beam path at an angle is the Swivel arm with scale. It enables you to swivel your Microbench construction in a controlled and measured manner. The swivel arm is screwed directly onto an optical table or fixed to a mounting plate using M6 or 1/4" screws.

This example demonstrates an assembly for measuring the angular intensity distribution of dispersion plates. The dispersion plate, backlit by an LED of the HighLED series, is fastened in a Microbench mounting plate in the center of rotation of the swivel arm. A photodetector is also installed in a mounting plate and mounted onto the mobile part of the swivel arm. This enables quick and easy measurement at any angle.

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What Adapters are Available?

Because there are so many different optic diameters and threads on the various construction elements, it is not always possible to install the desired components directly into the Microbench, Nanobench or Tube Mounting System C. Please take a look at the table below to find the adapters you need.

For details and technical drawings, refer to the listing for the specific system in the LINOS catalog or on our website www.qioptiq-shop.com.

| Artikel-Nr. | Außen-ø (mm) bzw. Gewinde | Inneanschluss bzw. ø (mm) | Länge (mm) |
|-------------|---|---|-------------|
| Order No | Outer ø (mm) or thread | Inner connector or ø (mm) | Length (mm) |
| G067023000 | 12.5 | FC | 8 |
| G067024000 | 12.5 | SMA | |
| G050310000 | 16 | FC | 8 |
| G050311000 | 16 | SMA | |
| G038855000 | 25 | 23.2 | 10 |
| G050305000 | 25 | 16 | 10 |
| G061626000 | 25 | 15 | 12 |
| G065042000 | 25 | 25 | 20 |
| G067054000 | 25 | FC | 10 |
| G067055000 | 25 | SMA | 10 |
| G067065000 | 25 | FC 3.5° | 10 |
| G038857000 | 30 | C-mount/FC 3.5° | 6.3 |
| G038851000 | 30 | C-Mount | 10 |
| G061602000 | 30 | 25 | 16 |
| G061634000 | 30 | 23.2 | 20 |
| G061658000 | 30 | C-Mount | 16 |
| G061637000 | 35 | 30 | 16 |
| G061674000 | 50 | For MB mounting plates | 13 |
| G061601000 | 25/30 | Solid material | 20 |
| G061621000 | 25/30 | W0.8"x1/36" | 20 |
| G061622000 | 25/30 | 16 | 20 |
| G061623000 | 25/30 | 19 | 20 |
| G061636000 | 25/30 | 23.2 | 20 |
| G061664000 | 25/30 | C-Mount | 22 |
| G061661000 | 25/C-Mount | 16 / W0.8"x1/36" | 34 |
| G061659000 | C-Mount | Connector for Microbench rods | 21 |
| G061667000 | C-Mount | M23.2x0.75 | 8 |
| G061668000 | C-Mount | W0.8"x1/36" | 8 |
| G038878000 | M17x0.75 | C-Mount | 11 |
| G061638000 | M35x0.5 | 30 | 16 |
| G061662000 | W0.8"x1/36" | C-Mount | 14 |
| G061085000 | 40x40, Adapter plate for FLS 40 carrier | Solid material | |
| G061005000 | 40x40, Adapter plate | Connector for Microbench and Nanobench rods | 10 |
| G050309000 | 25 | 16 | 10 |
| G060407000 | 30/C-Mount | M23.2x0.75 | 13 |
| G061168000 | 30 | 1.035"-40 UN | 10 |
| G173013000 | 30/50 | C-Mount | 14.5 |

Selection of adapters from the table

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How Do I Mount My Assembly, e.g. on Optical Tables?

To secure the Microbench system on an optical table or a breadboard, base or mounting plates are needed. These are made of aluminum and are available in various sizes and models.

☐ Options such as Mounting base 80x15 and Base plate 80x40, for example, let you install Microbench structures on table tops independent of existing hole patterns. Another option is Rod holder F, for mounting optical structures on surfaces with tapped-hole patterns. With these options, structures are mounted from the top with M6 screws.

- Fasten structures onto table tops independent of existing hole patterns with Mounting base 80x15 and Base plate 80x40
- Rod holder F lets you mount structures on surfaces with tapped-hole patterns

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Shortening Rod holder V, G06 1234 000

☐ Rod holder V can be individually shortened to install Microbench components at the beam height you require. These holders feature throughbore holes to mount them from above with M6 or 1/4" screws; the maximum beam height is 100 mm. To adapt the rod holder for your application, the unneeded length is sawed off the bottom of the holder.

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Advantages of Rod holder V:

- Additional adapter plates not required
- Stringent production tolerances ensure compatibility with other Microbench components; this is extremely difficult to achieve with parts you make yourself

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
☐ Small Microbench structures can be mounted on Mounting plate 220x40. This plate has countersunk bore holes for M6 screws in a 50 mm pattern, for mounting on table tops.

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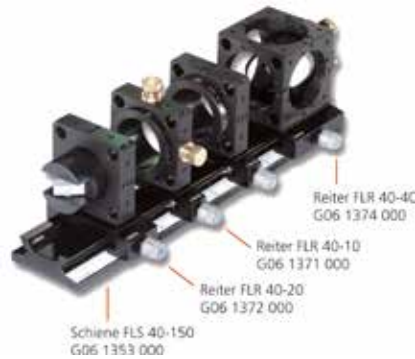☑ Microbench components can be mounted on the carriers of our Rail System FLS 40. This optical bench structure features high stability and, thanks to the slotted holes on the rail, can be fastened to tables and breadboards independent of existing hole patterns.

Flexible attachment options for optical tables with Rail System FLS 40

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☑ Nanobench structures can be fastened to optical tables or breadboards as well. Like the Microbench system, the mounting plates for Nanobench are made of aluminum and are available in various sizes and models.

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☑ The Tube system C also gives you several options for mounting components. Clamp holders 30 and 35 provide a wide range of mounting options. They can be mounted directly on an optical table or e.g. on the Base plate 80x40 with column.

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☑ The Clamping fork and Holder A14 serve to connect LINOS bench systems to optical tables. The slotted hole on the Clamping fork allows you to mount structures on optical tables without being restricted to a specific bore hole pattern. M6 or 1/4" threads are the only prerequisite.

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All components with M6 threads, in particular columns and support rods, are fastened to Holder A14 from below through a countersink. Simply slide the holder into the Clamping fork, position it as desired and then fix it in place from above through the slotted hole on the fork. The special fit of holder and clamping fork lets you lift them as a single unit to reposition the elements as desired. This is an exclusive design feature, called the easy lift system, and is especially practical when space is limited. In the example shown here, a support rod is clamped into Column 14, which gives you the advantage of an adjustable height.

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Tungsten-Halogen
Lamphouse
G06 5005 000

What Light Sources Are Available?

We offer a wide range of light sources, from bulbs to lasers, which are compatible with our bench systems.

☑ Use the Tungsten-Halogen Lamphouse for direct connection to the Microbench system. A concave mirror and an aspheric condenser are already included. On the front of the lamphouse a Microbench mounting plate ensures the direct connection to the Microbench system.

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LED Industry
Lamp IL100-c
G40 5817 000

☑ The LED Industry Lamp IL100-c is a versatile and robust LED Lamp for machine vision as well as for industrial applications. The eco-friendly and durable LED technology in combination with the optimized electronic allows for a safe and easy usage. The outer diameter of 25 mm enables the direct integration in Microbench mounting plates 25.

Assembly example: Microbench mounting plate fastened on a Flat Rail FLS 40 with Carrier FLR 40

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☑ The super-bright Microbench lights in the LED ML3 Series emit homogenous and flicker-free light. A digital potentiometer is used for processor-controlled adjustment of the brightness with no noticeable change in the color temperature. The LED ML3 Series Microbench Lights offer all the advantages of the latest in LED technology, as an energy efficient and environmentally friendly alternative to halogen lights.

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LED Lamp ML3 in Microbench mounting plate fastened on a Flat Rail 40 with carrier FLR

Several versions of white light are available in pure-white, natural-white and warm-white as well as UV-LED with 405 nm and 375 nm wavelength.

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The outer diameter of 30 mm fits the Microbench mounting plates. The combination with the Microbench and the Flat Rail System FLS enables perfect adjustability.

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☑ Another LED light source serie is our HighLED with adjustable brightness. The housing diameter is 30 mm for integration into Microbench components without additional adapters. Thanks to the integrated C-Mount connection threads, direct use in the Tube System is possible. We offer these LED light sources in four different wavelength ranges and with an extensive array of accessories.

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☐ Qioptiq offers helium-neon, or HeNe, lasers with Microbench-compatible outer diameters for the installation in the four-rod system. The laser beam is aligned and centered to fit the outer diameter of the housing. The lasers feature excellent directional stability, even when the ambient temperature is high. No X-Y adjustment is required following installation. Our HeNe-lasers have an outer diameter of 35 mm (direct installation) or 31.7 mm (using a Microbench adapter).

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■ With the Laser adapter C35, you can connect a HeNe-laser with a 31.7 mm diameter to the Tube system C.

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
The iFLEX-iRIS is a high performance solid-state laser system. All the smart control electronics and thermoelectric cooling (TEC) are integrated inside this small laser. The design includes automatic closed loop control in CW mode for excellent long-term power stability. Ultra-stable performance is also achieved when the laser is in Modulation mode. This allows you to achieve measurements that are precise and repeatable, day after day.

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
For more information, please refer to the iFLEX-iRIS brochure.

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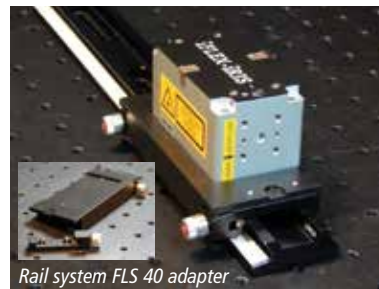✂ The iFLEX-iRIS laser modules have a bore-hole pattern which enables the mounting onto the FLS 40 rail system by using a FLS 40 adapter. Or you can directly fix the laser module onto an optical table by means of the iFLEX-iRIS Optical Table Mount.

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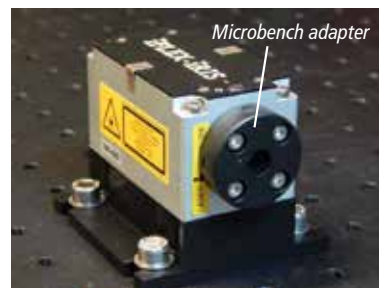☐ Standard Microbench mounting plates (such as G061010000) can be fixed to the Microbench adapter, which is screwed on the iFLEX-iRIS Laser module. Thereby, the iFLEX-iRIS lasers are compatible with Microbench setups. Of course, the optical axis of the exiting, collimated laser beam coincides with the central axis of the Microbench. This ensures ideal conditions for quick assembling with minimal adjustment effort.

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☒ A matching adapter also ensures the compatibility with the Nanobench system. The mounting plate N16 can be directly fixed to the Adapter.

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Please note that further Microbench compatibility is given with kineMATIX positioniers for kineFLEX optical fibers for fiber coupling iFLEX-iRIS Laser light. You will find more product information and further accessories for the iFLEX-iRIS solid state laser modules in our Q-Shop:

<https://www.qioptiq-shop.com/en/Lightsources-Laser/Laser/Diode-DPSS-laser/iFLEX-iRIS-accessories.html>

<https://www.qioptiq-shop.com/en/Lightsources-Laser/Laser/Diode-DPSS-laser/iFLEX-iRIS-accessories.html>

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How Can I Combine Systems?

☑ ⇨ With Optic Holder 63

With the Mounting ring for Microbench, the Microbench is compatible with Optic holder 63-M and can thus be integrated in this mount system as well.

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☑ ⇨ N

The Microbench adapter creates a connection between the Microbench and Nanobench systems by means of Mounting plate 30 and Mounting plate N 16. This adapter has through holes for Nanobench rods.

Another solution is the Microbench-Nanobench adapter. This adapter offers bore holes for Microbench as well as for Nanobench rods.

Adapter 16/25 enables the use of mounted Nanobench optics in Microbench mounting components with a 25 mm inner diameter.

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☑ ⇨ ☑

Adapters are available for the Tube system C as well, such as the Adapter C-Mount/Microbench, for example. You can combine the Tube system with the Microbench, lasers, CCD cameras, fibers and your choice of other instruments.

You want to construct or document the setup in your CAD software?

On www.qioptiq-shop.com we offer CAD files in dxf or step format. You will find the files for the respective products under the tab Docs + Drawings

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
Geräteflansch / Instrument flange

Mounting example for the instrument flange

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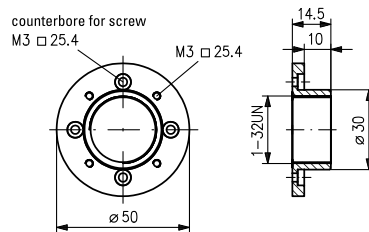
Beam deflection with Lees mirror mounts

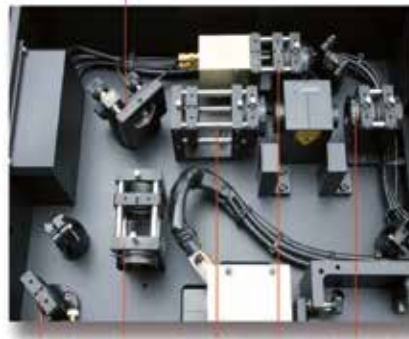Microbench components such as holders, mounting plates, rods, rotary supports, and so forth, in use

Beam deflector cube C-5; G06 5236 000

Clamp holder 30, G06 1243 000

How Do I Mount Bench Systems in My Instrument?

Our construction systems are ideal not only for laboratory structures on optical tables but they are also, thanks to their excellent stability, increasingly applied to measuring systems in both industrial environments and end-user instruments.

Components such as the Instrument flange or the Adapter FIC/Microbench allows you to connect modules made of Microbench or Tube System C components to your own instruments.

☐ As an example for the instrument flange, a laser beam coupling into a microscope is shown. The light of two lasers with different wavelengths is transmitted through fibers, collimated and consolidated through a beamsplitter cube with a structure made from Microbench and fiber optics components, and coupled in the microscope. Each collimator requires one adjustment option for superimposing the laser beams. The structure is connected to the microscope port using the Instrument flange.

☐ All Holders and Rod holders in the Microbench system have a centered M6 bore hole on the bottom and two M4 bore holes that are 20 mm apart. Rod holder F can also be mounted from above if desired, by means of an M6 through-bore hole. These bore holes make it easy to integrate Microbench structures in your own instruments. In addition the Microbench adapter gives you the option of connecting Nanobench components.

Our example shows the assembly of an Optical-Parametric Oscillator (OPO).

■ With Clamp holders from the Tube System, you can fasten any components with an outer diameter of 30 mm or 35 mm.

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SPCM with FC fiber adapter

SPCM with free space window

How do I connect Photon Counters to the Microbench?

Excelitas Technologies' SPCM-AQRH Single Photon Counting Module detects single photons over the wavelength range of 400 nm to 1060 nm with performance parameters superior to other solid state or vacuumtube based photon counters

The SPCM-AQRH uses a unique silicon avalanche photodiode (SLiK) with a circular active area, achieving peak photon detection efficiency greater than 70 % at 700 nm over a 180 µm diameter with unmatched uniformity over the full active area. A TTL level pulse is generated for each photon detected and the signal is available at the BNC connector at the rear of the module. The signal should be terminated into 50 Ω.

The photodiode is both thermoelectrically cooled and temperature controlled, ensuring stabilized performance despite ambient temperature changes. Operating temperature range has been increased and the module (case temperature) will function between 5 °C and 70 °C.

The module comes with a free space window, a FC fiber adapter, C-mount and tube cage mount bracket/adapter options.

Key Features

- Peak PDE: > 70% @ 700nm
- 180µm active area
- > 35 Mcps dynamic range
- Gated output
- Low after pulse probability
- High uniformity over large active area
- Unmatched linearity
- Free space window, FC fiber adapter, C-mount bracket, tube cage mount bracket options
- Timing resolution-optimized versions available

Applications

- LIDAR
- Quantum Cryptography
- Photon correlation spectroscopy
- Astronomical observation
- Optical range finding
- Adaptive optics
- Ultra-sensitive fluorescence
- Particle sizing
- Microscopy

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
spcm.europe@excelitas.com

SPCM-AQRH-WX-BR1 with bracket for 1" Tube (no threading) and 30mm Cage-System.
Mounting example with Microbench assembly

SPCM-AQRH-WX-BR2 with bracket for C-mount and 30 mm Cage-System:
Mounting example with Tube System C assembly

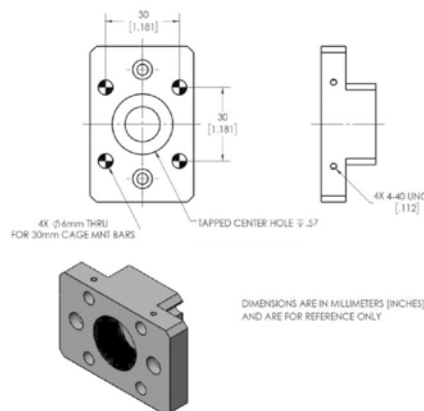Bracket BR1 adapter for Microbench Mounting plates 30 and Tubes with 1" Tube outer diameter

Optical Mount Bracket Options

There are two optical mounting bracket options available to the module. These mounting brackets are for industry standard C-mount and Cage Mount systems, such as the tube system C and the Microbench. In addition, the large portfolio of optical components can be integrated into these systems. Of course, other standard adapters (see page 27) can be attached to these SPCM adapters to ensure compatibility with a wide variety of applications.

Note: The mounting brackets are factory installed and cannot be removed by end users without risk of damaging the detector.

For more information, prices e.g. please contact:

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
How Do I Integrate Motorized Positioners?

Our series of high-precision piezoelectric positioning systems gives you a wide selection of miniature positioners in combination with the LINOS bench systems.

Positioner MS 30 is a miniature translation stage with an adjustment range of 8, 18 or 30 mm. The MS 30 is made of a warp-resistant aluminum alloy, has precision ball lining channels and can be mounted in XY or XYZ combination. The Miniature positioner MX 25 is a 3-axis adjustment unit with a 2 mm range.

The bore hole pattern of both positioners fit onto FLR 40 carriers (width: 34 mm and up). No further components are required for a combination with FLS 40 rails.

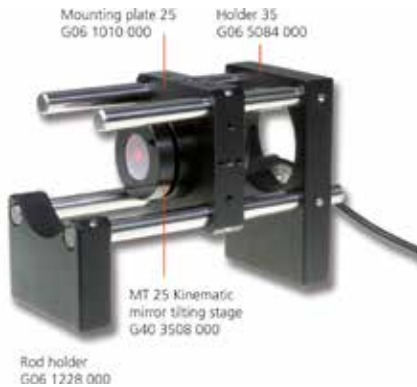
The Kinematic mirror tilting stage MT 25 with an outer diameter of 25 mm can be used in Microbench mounting plates and in cubes with the corresponding mounting diameter. Depending on your application, the MT 25 can be used for positioning mirrors or beamsplitter cubes.

A free through hole of 9 mm also permits the use of partially transmitting mirrors.

A manual controller is available. This controller CN 30 (G40 3513 000) enables the adjustment in three speeds using push buttons.

MT 25 Kinematic Mirror Tilting Stage in Microbench mounting plate

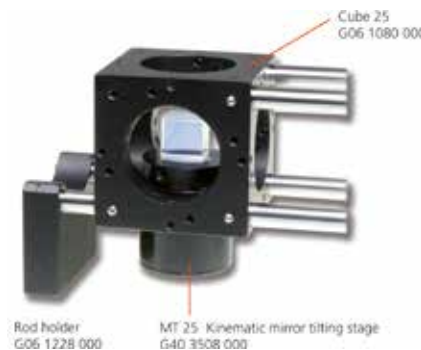Detailed information on all piezoelectric positioning systems is available in the Chapter entitled „Positioning Systems“ in the LINOS catalog and at www.qioptiq-shop.com.

Kinematic Mirror Adjustment MT 25 in Microbench Cubes

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Application areas:

- *Standard optics application in the lab*
- *Design allows applying a broad range of wavelength*

Q-Sets in our Shop:

Q-Sets – How to Save Time in the Lab

The Qioptiq Q-Sets provide you with complete, tested assemblies for your set-up.

A Q-Set is a set of components from the LINOS catalog which, when assembled, comprise a finished subassembly for typical optical setups in the lab. This saves your time right from the start, because you no longer have to search the catalog yourself for each individual component required for your optical or optomechanical system. All Q-Sets are thoroughly tested at Qioptiq. Thus we ensure the same high quality you get in all of our products. Each set is generally optimized for use with a wide range of wavelengths, to ensure broad compatibility.

You will find a variety of Q-Sets on the optomechanics homepage of our Q-Shop. All Q-Sets are available for immediate delivery.

Q-Sets are ...

- **Reliable:** performance tested in the Qioptiq labs
- **Compatible:** can be combined with established LINOS catalog components, especially with Microbench systems.
- **Time-saving:** complete assembly kit for standard lab applications
- **Beneficial:** less expensive due to set-price

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Q-Set Vertical Adjustment

This Q-Set is ideal for applications where the beam height needs to be adjusted in a range of 20-180 mm. Due to the small height of the lower beam-deflector cube the system can be easily adapted to the LINOS Microbench. This allows for the integration of a Microbench system in a configuration with a different beam height or vice versa. The maximum beam height is extendable via additional rods.

The mirrors are coated with a high quality silver layer with a reflectivity of more than 96% in the spectral range from 500 nm to 12 μm . They are suitable for applications with broadband light sources, as well as for femtosecond laser pulse applications. The surface-hardened post ensures a very high stability of the structure. The top mirror holder exhibits decoupled adjustment axes allowing a fast and precise alignment process.

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- Available for 488, 532, 633 or 780 nm
- Aperture: 9 mm

Q-Set Variable Attenuator / Beamsplitter mounted onto Mounting plate with a Mounting column 14 and Rod 14

Q-Set Variable Attenuator / Beamsplitter

The Q-Set fulfills two functions, one as variable attenuator where a rotation of the wave plate allows a continuous adjustment of the laser power in a system. Furthermore, it may be used as a variable beam splitter that splits the power between two orthogonal optical paths. The height of the Q-Set allows a combination with LINOS Microbench components. In that case the use of longer rods than the ones included is recommended. In addition, the mounting on posts as well as direct mounting onto optical tables is possible.

The zero order wave plates provide an increased wavelength, temperature and incident angle tolerance compared to wave plates of higher order. For the respective central wavelengths the wave plates come with an antireflective coating. The beam splitter cubes are equipped with a broadband antireflection coating.

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XY-translation stage:

- ± 1 mm adjustment range
- Decoupled adjustment axis

Z-Fine adjustment:

- Translation range: 1.5 mm
- Accuracy: 1 μ m
- Pinhole: 30 μ m diameter

Objective:

- 10x; NA 0.25; focal length 20 mm

Q-Set Spatial Filter

This Q-Set can be used to homogenize the intensity distribution of laser beams. Often dust and dirt particles on optical elements lead to aberrations. Also, the generation of the laser light itself may produce an unwanted intensity profile due to cavity properties. These aberrations can be reduced with the aid of the Q-Set Spatial Filter. The laser light is focused by an achromatic lens and transmitted through a pinhole which partially removes the inhomogenities. At the same time the Q-Set can be used to generate a point light source.

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- Attenuation up to 7 orders of magnitude
- Smallest increment in OD: 0.3

Q-Set Filter Cascade


The filter cascade provides a gradual attenuation of light over a huge wavelength and intensity range: the assembly of 5 neutral density filters enables an attenuation over more than 7 orders of magnitude for applications within the VIS and NIR wavelength range.

The carriers can be individually flipped in and out of the beam path and therefore allow the desired combination of filters for the respective application.

The Q-Set is equipped with neutral density filters with an OD of 0.3; 0.7; 1.3; 2.0 and 3.0. The assembly can be combined with other Microbench setups. Furthermore it can be directly mounted onto optical tables via the included FLS 40 rail.

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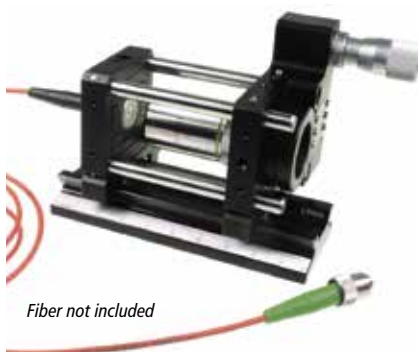Fle.X plate:

- Monolithic adjustment unit
- Decoupled adjustment axis
- Almost no thermal drifting: excellent long-term stability

Q-Set Fiber Coupler

The compact design of the Q-Set Fiber Coupler combines a flexible alignment with a high long-term stability. The z-Fine Adjustment unit holding the lens allows for a translation of 1 mm with 1 μm resolution and almost without play and tilting. The mounting and alignment of the FSMA fiber coupler is accomplished via a fle.X-plate which is a monolithic element with decoupled adjustment axes in the x and y direction. The respective adjusting range is ± 0.25 mm. Adjustment is carried out via tension and pressure screws, which are subsequently secured in position with lock-nuts for excellent long-term stability. The applicability for a broad wavelength range is enabled via the achromatic corrected lens.

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Q-Set Spectroscopy Liquids

The Spectroscopy Q-Set for liquids is designed for use with liquids in standard cuvettes with 12.5 mm edge length. Clamping the cuvette to the surface ensures reproducible measurement results.

The modular setup, with components from the Microbench and Flat Rail systems, makes it easy to adapt the structure to any application. The Q-Set includes a protective cap to prevent light scattering. The optic cell and fiber optic cables are not included in the set.

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Q-Set Spectroscopy Plates

The Spectroscopy Q-Set for plates makes it easy to perform transmission and reflection measurements on plates. Plates up to 10 mm thickness are clamped with protective plastic screws. The collimators can be adjusted to the desired height on the vertical rods, with a Fine Adjustment C30 integrated on one side.

The Q-Set consists of components of the Microbench and Flat Rail systems and enables individual adaptations and extensions. The foundation of the set is a Flat Rail 40-200 for easy attachment in devices for mounting on surfaces.

The fiber optic cables and the glass plate are not included.

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Q-Set Retroreflector

The Q-Set enables the rapid construction of a reflector. The angle between the mirrors can be varied almost arbitrarily and also the distance between input and output beam is variable, so that even smallest distances are adjustable. An angle fine alignment of the two deflection mirrors is possible via the adjustment screws. The mirror holders are mounted onto a vertically oriented flat rail via carriers. The rail in turn is mounted onto a clamping block attached to a solid steel column and therefore the structure provides an excellent stability. The distance between the mirrors can be adjusted along the rail to the desired position. The height of the entire system is varied via the clamping block.

In combination with the Q-Set Prism Stage, the Q-Set Retroreflector can be used to assemble a prism compressor.

Three M6x12 screws are required for assembling. Mirrors and the mounting plate are not included.

- Large variable deflection
- Height expandable by Column 25 (G061201000)
- Hardened column, diffusely reflective

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- Alignment with up to 5 degrees of freedom (3 rotation, 1-2 translation)
- High-precision positioning of prisms and beamsplitters
- For beam heights 12-21 cm (expandable)

Q-Set Prism Stage

The Q-Set allows for a highly accurate and flexible positioning of a prism for different beam heights. The prism table allows an adjustment of the prism in three degrees of freedom about a common pivot point. In addition, a horizontal alignment is possible via the supplied rail (coarse adjustment) and the linear stage (fine adjustment).

The mounting column with rod enables adjustment to beam heights 12-21 cm whereas the maximum height can be extended by using longer columns. An application example for this Q-Set is, together with another identical set or the Q-Set Retroreflector, a prism compressor.

The prism is not included in the delivery. For assembling the following screws are needed : Set screws M6x16 (1x), M2.3x5 (4x); Screw M6x20 (1x), M2.3x3 (4x) M3x6 (4x) M4x8 (4x).

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
The minimum focal spot radius w_f is yielded approximately by:

$$w_f \approx \frac{\lambda \cdot f}{\pi \cdot w_0}$$

In this equation, w_0 is the beam radius before a lens with focal length f .

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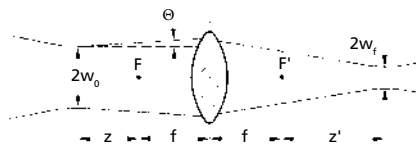Creating a Spatial Filter System

Scattered laser light can form interference-patterns with inhomogeneous intensity distribution caused by dust particles on optical surfaces. Such patterns can create problems with interferometric or holographic applications, because they form a confusing background to the desired information. Spatial filtering is a simple way to eliminate this interference and maintain homogenous beam intensity. The scattered light radiates outward from the laser beam in various directions and is thus spatially separate in the focal plane of a lens.

By centering a small aperture around the focal point of the direct laser beam, the scattered light can be blocked while permitting the direct beam to pass without interference. The result is a cone of light with an extremely homogenous intensity distribution that can be focused to a highly homogenous beam. As a compromise between ease of adjustment and complete spatial filtering, the best aperture diameter is approximately twice the $1/e^2$ diameter of the laser beam in focus, or roughly 1.33 times the 99% contour of the laser beam.

For many applications, particularly in holography, spatial intensity variations are unacceptable. Our spatial filter module is ideal for creating a clean Gaussian beam.

In this example, an Achromatic lens with $f = 40$ mm is mounted in the Z-Fine adjustment M and the laser beam is focused through the $50 \mu\text{m}$ pinhole. The Pinhole is in the Spatial filter module, which permits backlash-free XY adjustment thanks to a flexure hinge.


The focusing optics and pinhole size can be adapted to the particular optics used. The diameter of the pinhole should be approximately 30% larger than the calculated diameter of the focal point.

This application example shows one possible setup using Microbench components. A complete spatial filter module can also be found as a complete assembly kit (G06 0406 000) under the section Q-Sets (please look on page 38).

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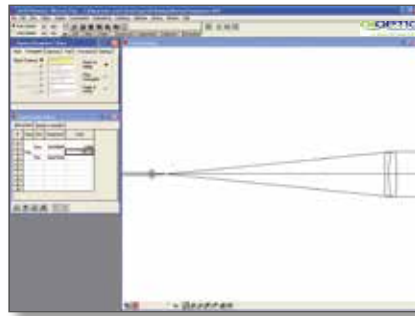Simulation of beam path with Qioptiq WinLens™

Creating an Expansion System

The Tube system C lets you set up simple expansion systems in accordance with the desired expansion ratio. With the Spatial filter module, you can add pinholes to the system as well.

Our example demonstrates a 16x expansion system. The structure of the optical system is simulated using the WinLens™ program from Qioptiq. The expansion system consists of two achromatic lenses, $f=10\text{ mm}$ and $f=160\text{ mm}$. These are positioned at a distance of their focal lengths, which forms a Kepler telescope. The pinhole is placed in the focal point.

For detailed information and a free version of the WinLens™ 3D software, point your browser to www.qioptiq-shop.com

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
This structure can also be realized with Microbench components. To do this, use the 4-rod structure, mounted Microbench lenses, a pinhole installed in a XY Adjustment plate, and the Z-Adjustment M.

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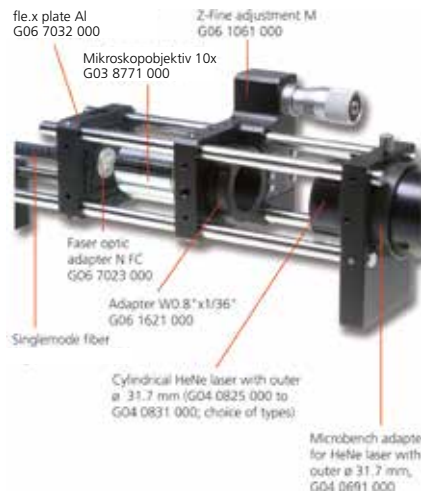You want to construct or document the setup in your CAD software?

On www.qioptiq-shop.com we offer CAD files in dxf or step format. You will find the files for the respective products under the tab Docs + Drawings

Creating a Fiber Coupler

There are a number of methods available for coupling light into a glass fiber, one of which is presented here. The task in this example involves fiber coupling a free laser beam into a singlemode fiber. The difficulty with singlemode fibers lies in their small core diameter which, depending on the wavelength used, may be just a few μm . The example illustrates how to solve the problem.

The fiber connector in the singlemode fiber is mounted in a special XY Adjustment plate adapter that permits fine adjustment in the range of $\pm 0.25\text{ mm}$ and provides excellent long-term stability. The XY Adjustment plate is a solid plate with structures generated by wire-electro discharge machining. The result is a broad solid-state spring joint that enables extremely precise XY translation with no tilting. Translation is achieved using tension and thrust screws secured against one another by lock-nuts to ensure long-term stability.

A Microscope lens is used for fiber coupling of the laser light. Proper adjustment of the lens along the Z-axis is indispensable for optimum coupling efficiency. In our example, we use Z-Fine adjustment M. With this high-precision, backlash-free positioner, the optimum distance from the microscope lens to the optical fiber is adjusted quickly, and high coupling efficiency is attained.

| Artikel-Nr. | Produkt |
|--------------|---|
| G06 7032 000 | fle.x Plate AI |
| G06 7023 000 | Faseradapter N FC |
| G03 8771 000 | Mikroskopobjektiv 10x |
| G06 1061 000 | Z-Feinverstellung M |
| G06 1621 000 | Reduzierring W08"x1/32" |
| G04 0691 000 | Mikrobänk Adapter für HeNe-Laser, Außen-ø 31.7 mm |

| Order No | Product |
|--------------|--|
| G06 7032 000 | fle.x Plate AI |
| G06 7023 000 | Fiber Coupler N FC |
| G03 8722 000 | Microscope Lens 10x |
| G06 1061 000 | Z-Fine Adjustment M |
| G06 1621 000 | Adapter W08"x1/32" |
| G04 0691 000 | Microbench Adapter for HeNe-Laser, outer ø 31.7 mm |

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A HeNe laser for this application example needs a laser head with an outer diameter of 31.7mm. Suitable HeNe series with further information can be found in the Q-Shop under:

Another, fully configured fiber coupling unit from Microbench components can be found in the Q-Sets section (G060408000 on page 39)

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
Setup of a Modular Flow Cytometer

How to rapidly prototype a complex photonic system by using the LINOS Microbench: A joint project of Excelitas, Qioptiq and Kinetic River.

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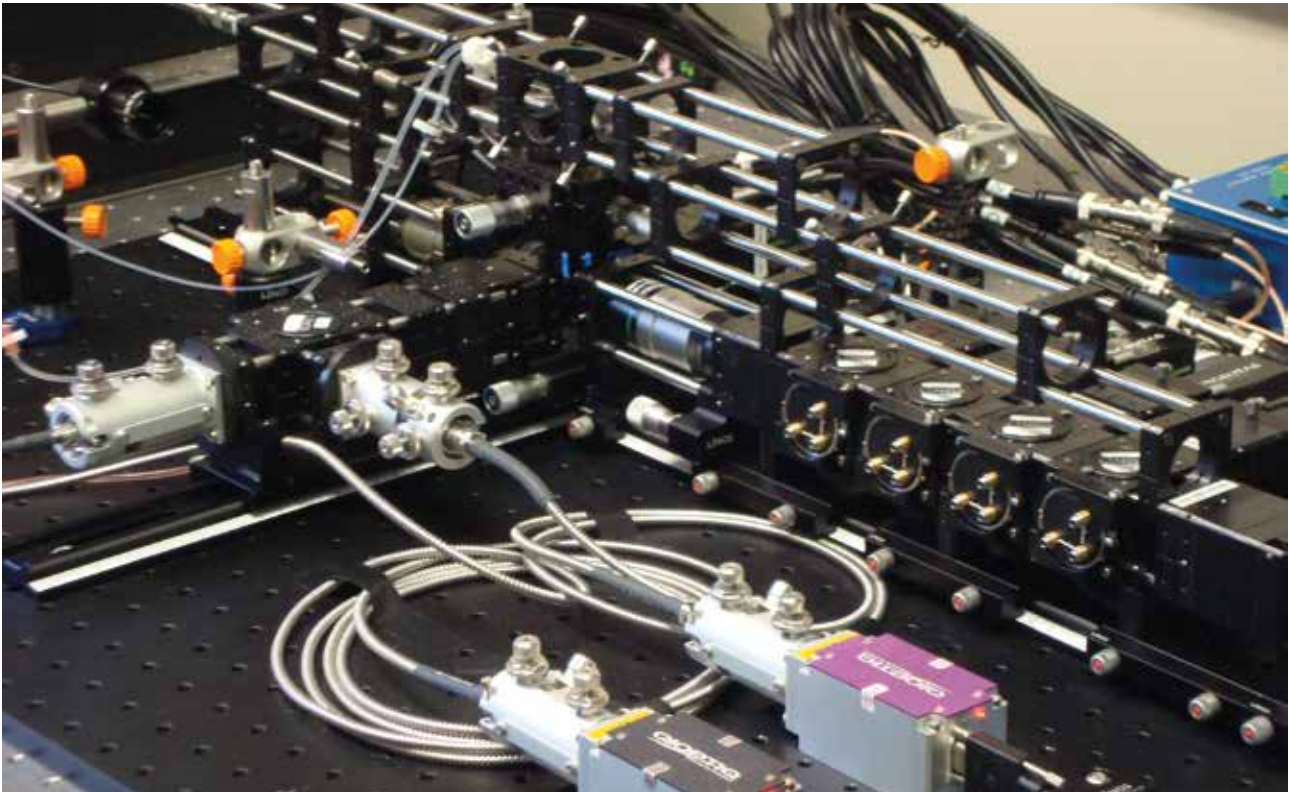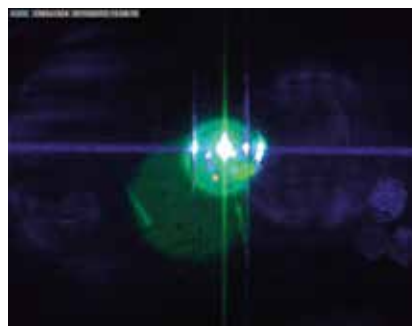
The violet laser light generates a fluorescent sample stream.


About flow cytometers

Flow cytometers are instruments that measure cells in a liquid suspension, characterizing cell populations using light scatter and fluorescence parameters. These life science instruments play an important research role in immunology and cell biology, and they are used in the clinical lab for diagnosis and monitoring of diseases such as HIV / AIDS and blood cancer.

In a flow cytometer, samples of cells, often mixed with antibodies and fluorescent dyes, travel through a fluidic path to a glass flow cell, where each individual cell is intercepted by light, generally from laser sources. Collection optics and filters direct the scattered and fluorescent light to the detectors. The signals are digitized for visualization and analysis.

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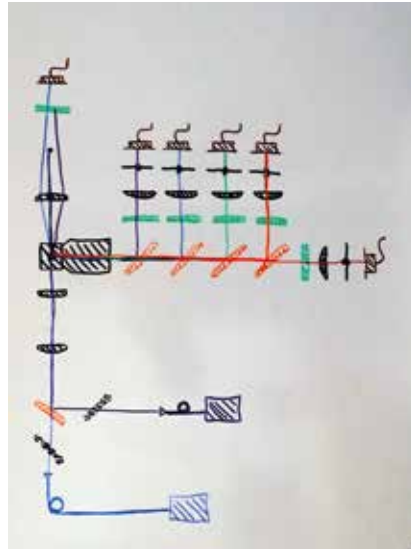
Kinetic River Corp., a leader in biophotonic design, consulting, and product development, designed and built this 4-color flow cytometer prototype from scratch in less than 8 weeks, using almost exclusively off-the-shelf Excelitas and Qioptiq photonic components. The development of this fully functional instrument prototype demonstrates the ease and expediency of designing and building life science and analytical instruments from the ground up using the LINOS Microbench system.

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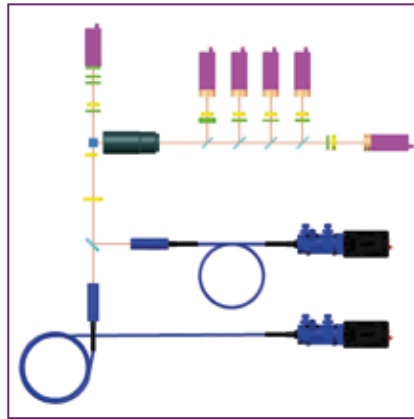Arrangement of the required Qioptiq iFLEX-iRIS solid-state lasers (including fiber couplings) for the light scattering at the sample particles and for the excitation of the dyes coupled thereto. Six Excelitas LynX silicon photomultiplier modules enable the detection of forward scattering and sideways light scattering as well as the realization of 4-fluorescence channels. The entire optical system also includes various color filters, beamforming lenses, collecting lenses, diaphragms, dichroic beamsplitters, and a microscope objective.

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
Components

- 2 fiber-coupled Qioptiq iFLEX-iRIS™ excitation lasers
- 6 Excelitas LynX™ SiPM detection modules
- Qioptiq kineMATIX™ 4-axis micromanipulators
- Qioptiq OPTEM microscope objective
- LINOS Microbench and rail systems

Applications

- OEM instrumentation development
- Platform design
- Rapid prototyping for feasibility testing

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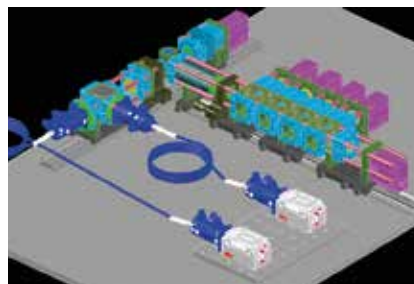
Lasers, detectors, and optics integrated in Microbench components and mounted on the Railsystem FLS 40.

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
Developed framework for mounting the flow cell (1), for adjusting the optical path of the laser light (2) and the scattered light (3) and for attaching the monitoring digital microscope (4).

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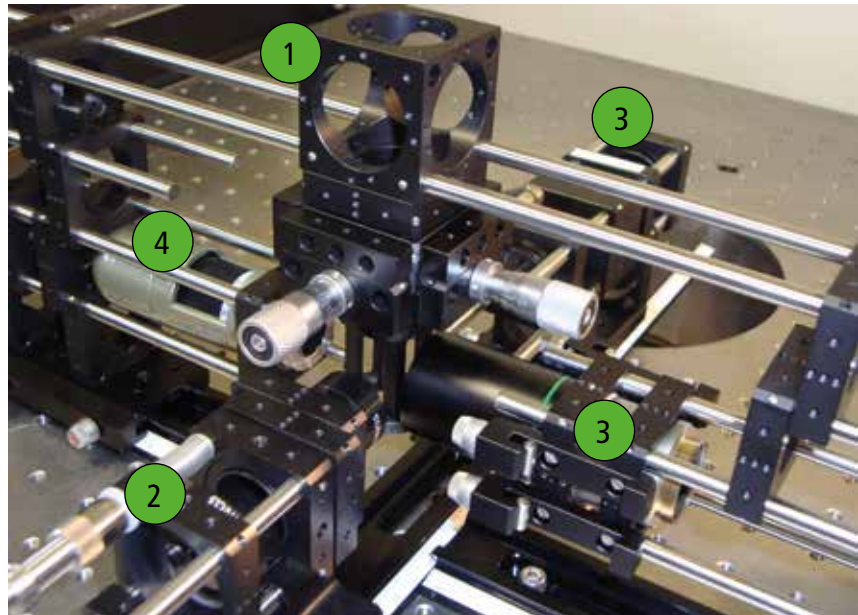*8 weeks later:
Fully assembled and tested prototype presented later at Photonics West, California, in February 2015.*

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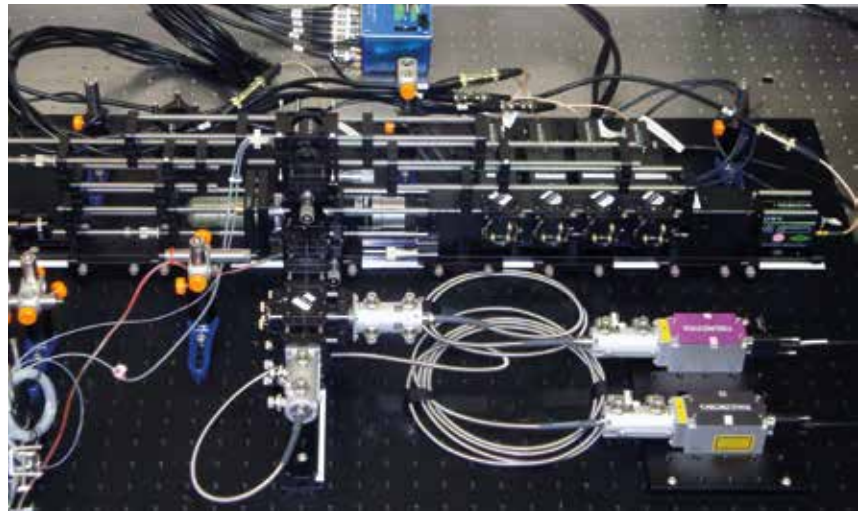Conclusion

The open, modular Microbench architecture of this prototype enables great flexibility in system layout and component selection for designers. Using off-the-shelf components significantly reduces the development time and costs.

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Have a look at the YouTube video for more information about the design, modularity and the functionality of the finalized prototype.

Discover the Q!

Qioptiq stands for products of the highest quality, whether off-the-shelf optical and mechanical equipment or customerspecific components and modules. Qioptiq designs and manufactures photonic products and solutions that serve a wide range of markets and applications in the areas of medical and life sciences, industrial manufacturing, defense, and aerospace, as well as in research and development.

Since 2013, Qioptiq has been an Excelitas Technologies company. Excelitas supplies innovative solutions for illumination, detection and other high-performance electrooptic technologies. The portfolios of these two companies complement one another perfectly, and together they employ more than 6000 people worldwide.

Photonics for Innovation

Contact Qioptiq today:

North America
+1 (800) 429 0257

Asia/Pacific
+65 64 99 77 77

Europe
+49 (0) 551 69 35-0

qioptiq.com
qioptiq-shop.com

© 2018 Qioptiq Photonics GmbH&Co KG. All rights reserved.

Qioptiq reserves the right to change this document at any time without notice and disclaims liability for editorial, pictorial or typographical errors.

O-CA_BR-Microbench-DE-EN_2018.02